


BED SOLUTIONS ICU BEDS

SEMI MOTORIZED
MANUALLY OPERATED

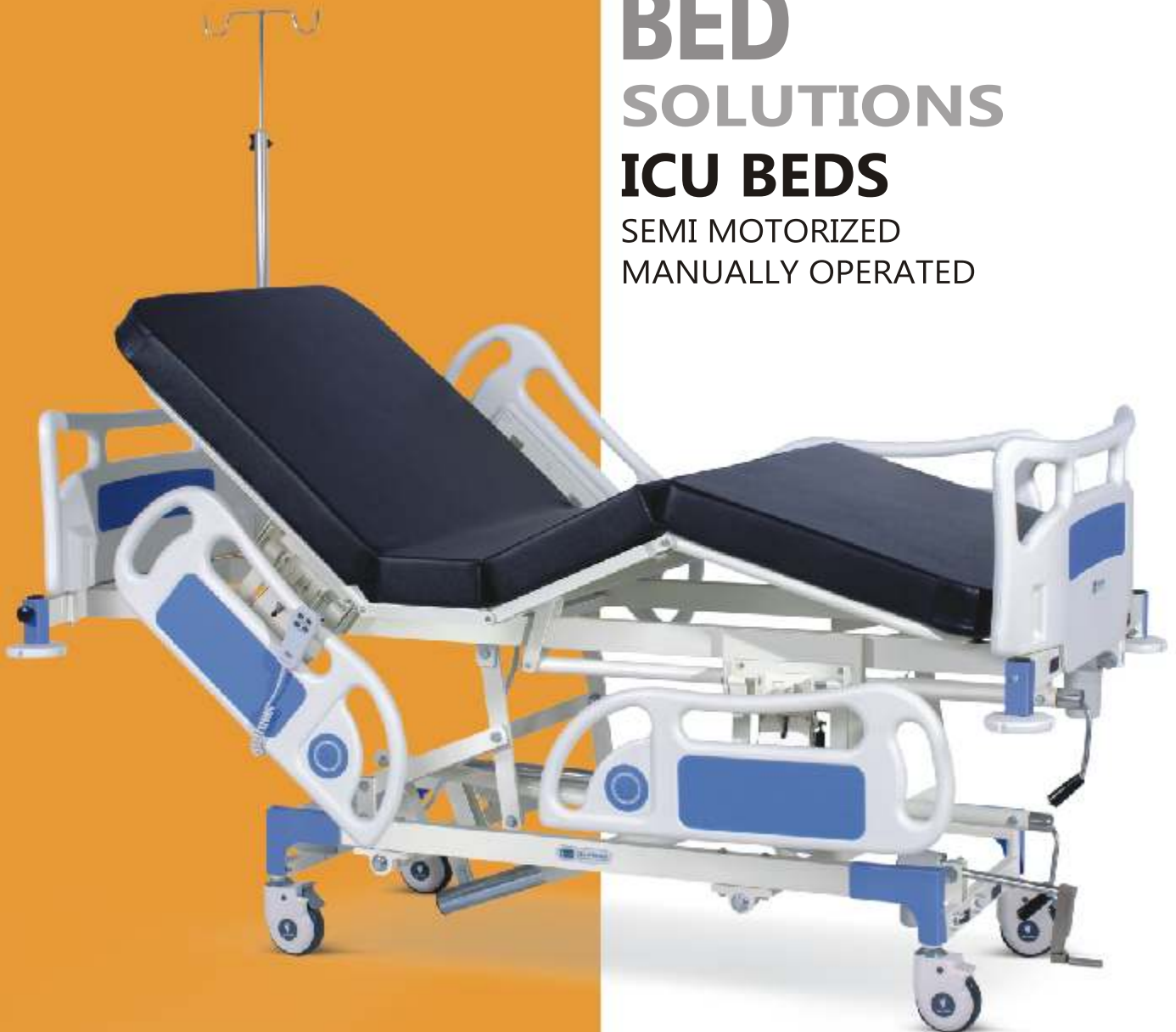

Mi-9020 AX

**Semi Motorized ICU Bed
with Backrest and
Height on Motors**

Introduction:

Medimek is a trusted brand in the medical fraternity over a decade. Our singular and simple philosophy is to cater to the needs of all our customers with a comprehensive range of excellent quality, user friendly hospital furniture and allied products, offered in a broad spectrum of exciting shades, aesthetic patterns and creative designs and backed by prompt after sales services.

Beds are universal and essential for ensuring successful medical treatment to the patient in any hospital and medical institution. Adding option of motorized function enhances patient comfort and reduces stress of the care givers.

ICU Bed is a holistic product designed with excellent features for patient positioning together with mobility, comfort and luxury. A composite hospital bed, it is acclaimed for Emergency care and Critical care because of its ability to offer all primarily important positions.



Mi-9020 BX

Semi Motorized ICU Bed with Backrest and Height on Motors
with Polymer head and foot boards and Collapsible railings



Mi-8001 AX

Manually Operated ICU Bed
with Polymer head and foot boards and Polymer railings



Mi-8001 BX

Manually Operated ICU Bed
with Polymer head and
foot boards and
Collapsible railings



Mi-8002 CX

Manually Operated ICU Bed
with SS head and
foot boards with Metal panels
and Collapsible railings



Backrest



Upper
leg



Height



Trende-
lenburg



Reverse
Trende-
lenburg



Fixed SS
handle



Dual
function
detachable
handle

Options Available :

Semi Motorized ICU Beds:

Electrical and Manual adjustments:

- Only backrest electrically adjusted
- Backrest and height / upper leg electrically adjusted OR
- Backrest, upper leg and height electrically adjusted with user friendly handset
- Trendelenburg and reverse trendelenburg and any non-motorized adjustments on separate screw mechanisms provided with handles
- Electrical adjustments available either on actuator with manual override system OR actuator with battery back-up facility

Height adjustment:

- Height adjustment (approx. without mattress): 515 mm to 720 mm

Manually Operated ICU Beds:

Manual adjustments:

- Backrest, upper leg, height and trendelenburg and reverse trendelenburg adjustments on separate screw mechanisms provided with handles

Height adjustment:

- Height adjustment (approx. without mattress): 540 mm to 740 mm

Common Features:

- Four section uniformly perforated mild steel CRCA sheet top
- Top frame and base frame made from mild steel tubes
- Telescopic stainless steel IV pole with four locations
- Four heavy duty castors 125 mm dia, two with brake
- Rubber buffers at four corners
- Pretreated and powder coated finish
- Supplied in SKD condition

Head and Foot boards :

- Polymer head and foot boards / Stainless steel head and foot boards with Metal panels

Bed dimensions :

Option 1:

- Bed frame approx. dimension : 2080 mm L x 915 mm W [82"L x 36"W]
- Overall approx. dimension (buffer-to-buffer) : 2180 mm L x 1015 mm W

Option 2:

- Bed frame approx. dimension : 1980 mm L x 840 mm W [78"L x 33"W]
- Overall approx. dimension (buffer-to-buffer) : 2080 mm L x 940 mm W

Bed frame dimension varies with Polymer railings

Optional Accessories:

Mi-1002

Polymer railings with degree indicator (set of four)

Mi-1004

Polymer railings without degree indicator (set of four)

Mi-1007

Collapsible railings with SS tubes and MS brackets (pair)

Mi-1021

Four section 100 mm thick PU foam mattress covered with rexine

Mi-1026

Pillow covered with rexine

Mi-1066

Heavy duty Telescopic Infusion Pump pole - SS

Mi-1076

X-ray permeable backrest with SS cassette holder

Mi-1077

Lifting Pole - SS

Mi-1082

Chart Holder

- Accessories, colors and designs shown in the brochure illustrations may not be the part of standard product configuration. • We reserve the right to make changes to designs, specifications and illustrations without prior notice. • Reproduction of this brochure in any form is liable for legal action.